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
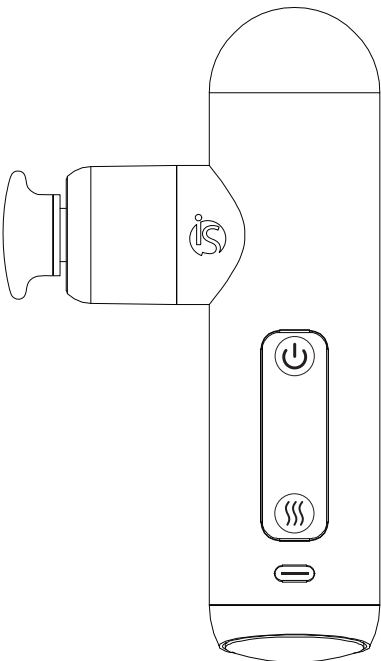
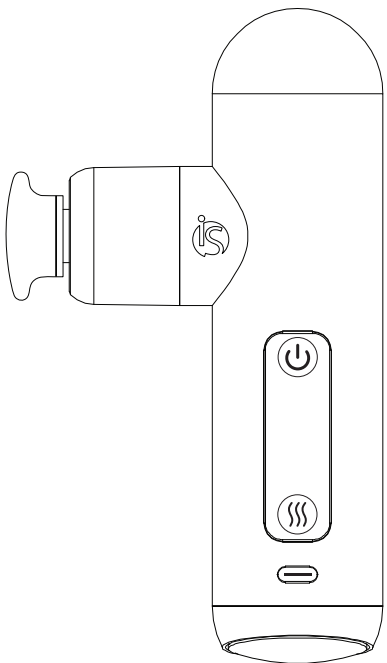
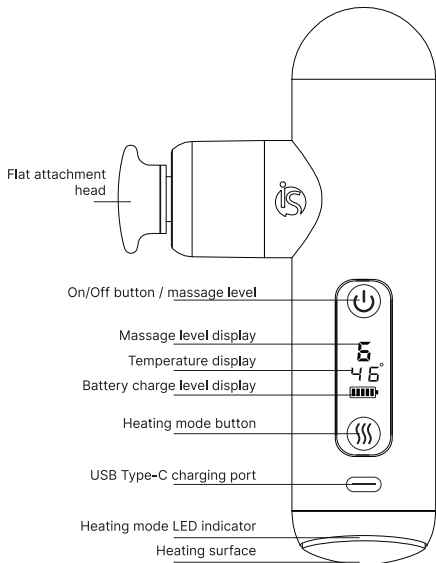
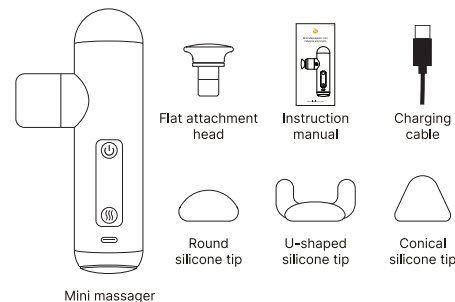
Using the Product


Assembly

With the massager switched off, insert the flat attachment head.

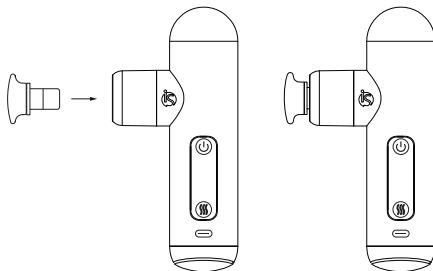Attaching the silicone massage tips

With the massager switched off, insert the desired silicone massage tip into the flat attachment head.

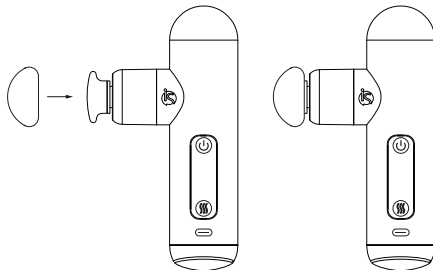Flat head: suitable for relaxation, especially the abdominal area, pectoral muscles, biceps, and femoral triceps.

Round head: suitable for large muscle groups, including arms, waist, back, buttocks, and thighs.

Conical head: suitable for small areas such as the soles of the feet, palms of the hands, the Achilles tendon, and the rhomboid muscles.

U-shaped head: suitable for trapezius muscles, areas around the Achilles tendon, ankle muscles, and the lower part of the leg.

Using the Product

On/Off

Hold the On/Off button  for 3 seconds to switch the massager on or off.

The massager will automatically switch off after 15 minutes of operation. If you wish to continue using it, switch it on again.

Massage Mode and Intensity Levels

The massager will start operating in massage mode at the first intensity level, visible on the screen. Six levels are available, which you can select sequentially by pressing the  button.

	1700 RPM		2300 RPM
	1900 RPM		2500 RPM
	2100 RPM		2800 RPM

Heating Mode

The massager has a surface with heating capability, located in its lower area.

To activate the Heating function (after switching on the massager), hold the Heating Mode button  for 3 seconds.

Heating will start at 38 °C, shown on the temperature display. The LED indicator will light up green. You can then switch sequentially to other heating levels (42 °C and 46 °C) by pressing the  button again. The LED indicator will turn blue at 42 °C and red at 46 °C.

To switch off, hold the  button for 3 seconds. (The massager will turn off.)

Note: Massage mode and Heating mode cannot operate simultaneously. When Heating mode is switched on, Massage mode switches off.

To return to Massage mode, hold the On/Off button  until the massager switches on again at intensity level 1.

Using the Product

Battery and Charging

You can check the battery charge level on the screen. The indicator has 5 bars , each corresponding to 20% of charge.

When the charge level is low, the indicator will show a single bar flashing .

To charge the battery, switch off the massager and connect the included charging cable to a power source (5V/2A max) and to the Type-C charging port on the massager.

Do not switch on the massager while charging.

During charging, the bar corresponding to the current level will flash .

When all bars are lit, charging is complete .

Charge the massager before first use, or if you haven't used it in 12 months.

Specifications

Charging time	3 hours
Battery capacity and type	1800 mAh, lithium-ion
Charging interface	USB Type-C
Autonomy	3–5 hours
Power	13.3 W
Dimensions	90×40×155 mm
Weight	400 g
Noise level	45–50 dB
Materials	Aluminium alloy
Operating temperature	10–40 °C

Safety Instructions

- Always use a certified 5V charger/adaptor in good condition
- Do not use fast-charging adapters
- Charge the product at ambient temperatures between 10 - 40 °C
- Do not heat the device above 70 °C
- Avoid proximity to fire, water, and direct sunlight
- Switch off the product after use
- Do not use while charging
- Keep out of reach of children
- Do not disassemble or attempt to repair the equipment
- In case of damage, stop using immediately
- Do not use on injured or bleeding areas of the body
- Do not use on areas of the body with recent scars
- Do not use if you have a pacemaker/defibrillator implanted
- Do not use on the head, on bony prominences, or on the abdomen during pregnancy
- Do not use on areas of inflamed skin or skin with diseases
- Do not use in cases of circulatory diseases or cancer
- Do not use in cases of acute illnesses (fractures, ruptures, or others)
- Do not use on vital organs or in the genital area
- Do not use on joints or bones
- Do not bring hair close to the massager to avoid it getting caught and causing injury
- During use, do not touch the back of the attachment head to avoid injury
- The use of the massager does not replace or complement treatments recommended by healthcare professionals
- In case of doubt regarding the use of the device, consult a healthcare professional

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